

Assistive Listening Devices

Alarm Clocks

Telephones

Television Listeners

Personal Listening Devices

Loop Systems



Ampli-DECT 350 Cordless Phone

The Ampli-DECT 350 was developed specifically to be used with hearing aids in the 'T' position. Features up to 30dB amplification, tone control, caller ID, 4 memories and last 10 outgoing calls redial feature.



ACC390B Ampli-DECT 350 Cordless Phone

Wake 'n' Shake Alarm Clock



This stylish, large LED display digital alarm clock wakes even the heaviest of sleepers using sound, strobe light, and a vibrating pad (included). Combination options include noise & flash, noise & vibrate or flash & vibrate. Mains powered (UK plug, supplied).



ACC369A Wake 'n' Shake Alarm Clock

CL-10 Corded Telephone



Telephone for the hearing impaired, featuring 2 levels of hearing amplification; standard, medium (+10dB), and high (+20dB), as well as 3 speed-dial numbers. It also features a hearing aid compatible inductive coupler, ringer on/off feature and indicator light for incoming calls. The CL-10 is also wall mountable.



ACC365 CL-10 Corded Telephone

CL-1 Telephone Ring Flasher



Provides a light and sound indication for incoming telephone calls. Mains powered (supplied), plugs directly into a telephone jack and a mains power socket. Features sound and strobe flash options, as well as volume and tone controls. The CL-1 is wall mountable and includes a 2-way phone socket adaptor.



ACC367 CL-1 Telephone Ring Flasher

CALL20 Tel/Door Ring Flasher



Never miss a door bell or a telephone call again with the AmpliCALL20. When your door bell/telephone rings, this doorbell/telephone ring indicator will ring very loudly and flash to let you know there is someone at the door/or that there's an incoming call. Features extra loud ringer, adjustable volume, tone control and 4 alarm options. Wall or desk mountable.



ACC367C CALL20 Tel/Door Ring Flasher

Assistive Listening Devices

Mr Rooster Alarm Clock

Alarm clock that generates a very loud rooster call and features a vibrating pad for use under the pillow. Designed with the hearing impaired in mind. Supplied with 4 x AA batteries.



ACC212B Mr Rooster Alarm Clock

Audioport A200



The A200 is a stereo personal sound amplifier for people who may sometimes need to rely on assistance for their hearing. Due to the stereo amplification via two integrated electret microphones, the user can clearly determine the direction of the sound source within an 8' to 10' radius. Dynamic headphone transducers ensure a natural and transparent sound reproduction with optimum speech intelligibility.



ACC371 Audioport A200

ACC371B Replacement Battery

IQ Buds BOOST

IQ Buds BOOST is the most advanced assistive listening device on the market today. Built with all the great features as our award winning IQbuds; we've added Ear ID™, a clinically backed audiometric hearing assessment that calibrates the IQbuds to your unique hearing profile.

Features

- Ear ID™ - Self Fit, Self Assess, Auto Calibrate
- SINCTM™ - Blend and Augment ambient noise
- High Fidelity Audio
- Hands-free calling
- Tap-touch control with Siri and Google
- Customizable tips for a secure fit; Includes Comply foam tips
- Best-in-class battery life with 20 hours on-the-go audio streaming, and 32 hours assistive audio



IQBUDSBOOST Intelligent Wireless Hearing Buds



New Product



Sennheiser® TV Listener



The Sennheiser® 830-S is an infra-red stereo TV listening systems. It enables listening anywhere in the room without the need for cables and without disturbing others.

Directly connectable to a TV, Hi-Fi system or radio, it guarantees outstanding sound quality and optimised speech intelligibility, thanks to switchable compression and treble emphasis. Ergonomic controls and visual signals for the individual control processes means that operation couldn't be easier.

Available in inductive neck loop style.



ACC416A Sennheiser® Set 830-S

ACC371B Replacement Battery For 830 Models

Audimed Communicator 2

The Audimed Communicator 2 is a versatile, lightweight listening device that will provide amplification to help those with a moderate hearing loss, who do not use a hearing aid.

Typical Hearing Situations

Allows the hearing impaired to communicate easily without using a traditional hearing system.

Listen to television without disturbing others.

Great for boosting sound of speech in meetings or at the theatre.

Ideal for short term use where the user may not have access to their hearing system (e.g. hospital)

Features

Combined ON/OFF and TONE switch

(Off | Normal Response | High Tone Response)

Built-in microphone

Volume control

Lightweight detachable stethoclip headband with universal fit earpieces

Available in 5 colours.

Specifications

Battery life: 250-500 hours



Battery type: Size 675 zinc air

Output: 128dB

Gain position 1: 56dB (wide band)

Gain position 2: 45dB (low frequency cut)

Frequency Range: 200-6500 Hz

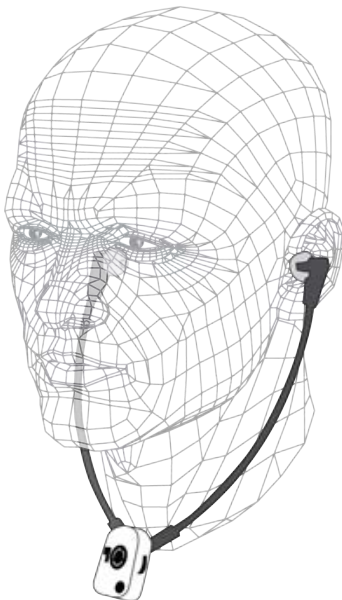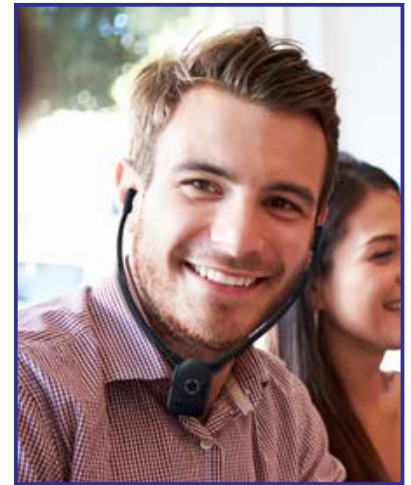

COMM-2-BK	Black
COMM-2-BL	Blue
COMM-2-CL	Clear
COMM-2-RD	Red
COMM-2-WH	White

Inductive Loop Listener

The Loop Listener is a compact listening device for interaction with hearing loop systems, It enables people to listen to a hearing loop system as if they were wearing a hearing aid switched to the "T" position. It is ideal for testing loop systems.

Applications


- To determine if a hearing loop system is functioning
- For someone with hearing difficulty who does not wear a hearing aid
- If a hearing aid does not have a "T" function

Features

- Flat frequency response
- Low cut filter
- More than 100 hours life when used with 2 x AAA batteries

What's Included

- The Loop Listener
- Headphones
- 2 x AAA batteries

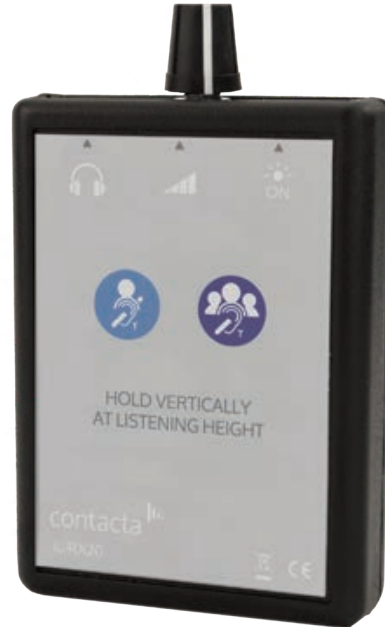

ACC338E Loop Listener

Portable Loop System

The Portable Loop is a versatile solution, providing mobile desktop assistance for your customers. It can be moved between rooms as required and is used for conversations between two people. The single control unit is housed in robust ABS plastic and incorporates the microphone, amplifier, hearing loop aerial and rechargeable battery system.

Applications

- Reception areas
- Shops
- Banks
- Information centres

Features

- Automatic gain control
- Compression & noise gate
- Constant-current loop drive
- LED indicator shows loop drive & battery status
- Magnetic field strength distance approx 1.2m
- Can be used with boundary microphone

What's Included

- The Portable Loop
- Power supply
- IEC lead
- Hearing loop signs
- Portable loop shelf



ACC338A Portable Loop System

ACC338B External Microphone

ACC338L UK Power Supply

Under The Counter Loop Kit



New Product

The “Under the Counter” loop kit provides clear communication with hearing aid wearers, especially where high background noise is a problem or clear speech is restricted. It has a flexible aerial that can be bent to shape, allowing it to be installed into many different positions.

The kit includes our halo indicator noise cancelling microphone as standard for high performance and self-testing. A kit containing a small omni-directional mouse microphone is also available.

Applications

For use in any situation where one to one communication takes place at a counter, booth or desk, such as in banks, post offices, supermarkets, theatres, airports and information centres.

Features

- Microprocessor control
- Energy-saving “green” standby mode is achieving a power saving of up to 80% during quiet periods
- Advanced audio processing features AGC, compression, and noise gate
- Constant-current loop drive
- Continuous self-testing
- Isolated audio line input allows the driver to take alternative audio sources (optional)
- Power and LED status indication



ACC338N Under The Counter Loop Kit - Halo Mic

ACC338P Under The Counter Loop Kit - Mouse Mic

Home Loop Kit

The “Home Loop” is a sleek and powerful hearing loop driver that connects to televisions and other audio devices. It enables hearing aid wearers to listen wirelessly directly from the source, providing the best possible sound quality. Integrated clear audio processing gives the user outstanding sound quality, enhancing speech and music clarity.

Applications

Suitable for home and business use, in small rooms where a television or sound system is present such as a living room.

Features

- Quick setup & simple connection to either a room loop, or a loop pad underneath a sofa or chair
- Optical IN for clearer sound less prone to interference
- Adaptable to suit different room layouts
- Connects to multiple devices
- Constant current for consistent audio performance
- Separate volume control
- HF compensation



ACC338M Home Loop Kit